


# Orbitrap IAPI Online licensing guidance material

 The world leader in serving science

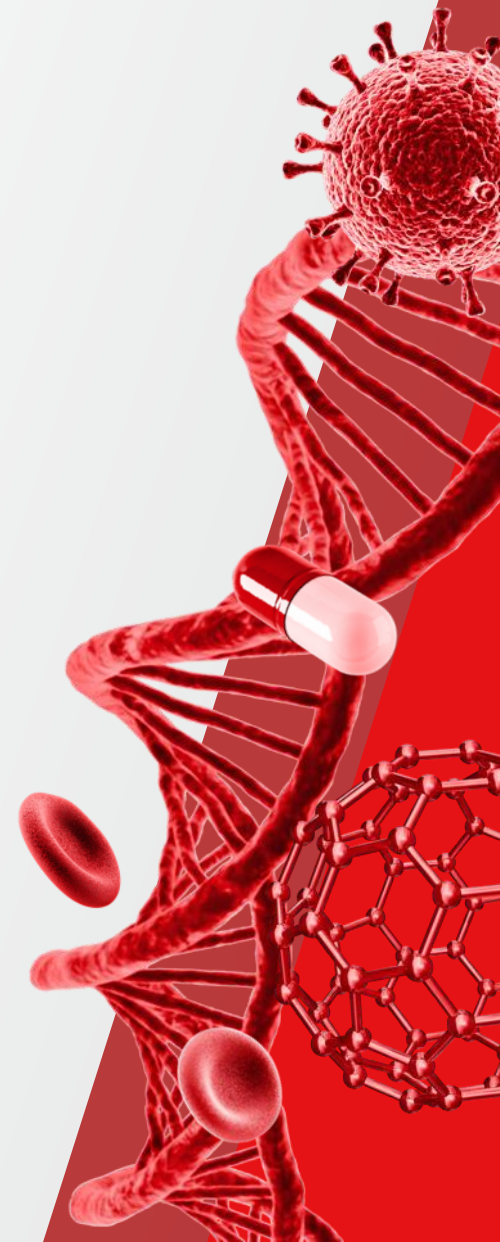

**1. Connecting system to Almanac**

**2. Online IAPI licensing**

# Connecting your system to Almanac

Learn more about Almanac at [www.thermofisher.com/almanac](https://www.thermofisher.com/almanac)

# Connect to Almanac via [www.thermofisher.com](http://www.thermofisher.com)

## Create an Account

1. Go to [www.thermofisher.com](http://www.thermofisher.com)
2. Click on 'Sign In' to the top right corner
3. Select Register to create an account

← → ↻ 🔒 [thermofisher.com/de/en/home.html](http://thermofisher.com/de/en/home.html)

1

ThermoFisher  
SCIENTIFIC

Popular Applications & Techniques Shop All Products Services Support

Search All

Search by catalog number, product name, keyword, application

3

ThermoFisher  
SCIENTIFIC

### Why you should create an account

Quickly and easily register to take advantage of these benefits:

- Obtain account-specific pricing and online quotes
- View and track existing or past orders and quickly reorder
- Join our loyalty program to receive a free, full-size sample
- Collaborate via a shared shopping list
- Shop the online scientific Services Marketplace
- Utilize 1TB of free data storage, scientific analysis apps, and peer collaboration tools

First name: \*

Last name: \*

Email: \*

This will be your username

Password: \*

Password strength:

Sign In ▾ Quick Order 

Sign In

Register

Account

Check Order Status

Aspire Member Program

Connect: Lab, Data, Apps

Custom Products & Projects

Instrument Management



Services Marketplace

2

# Connect to Almanac via Xcalibur

## Register the Almanac Agent

1. Open the Xcalibur Roadmap View
2. Click on the 'Cloud' icon to access Almanac
3. Launch registration wizard

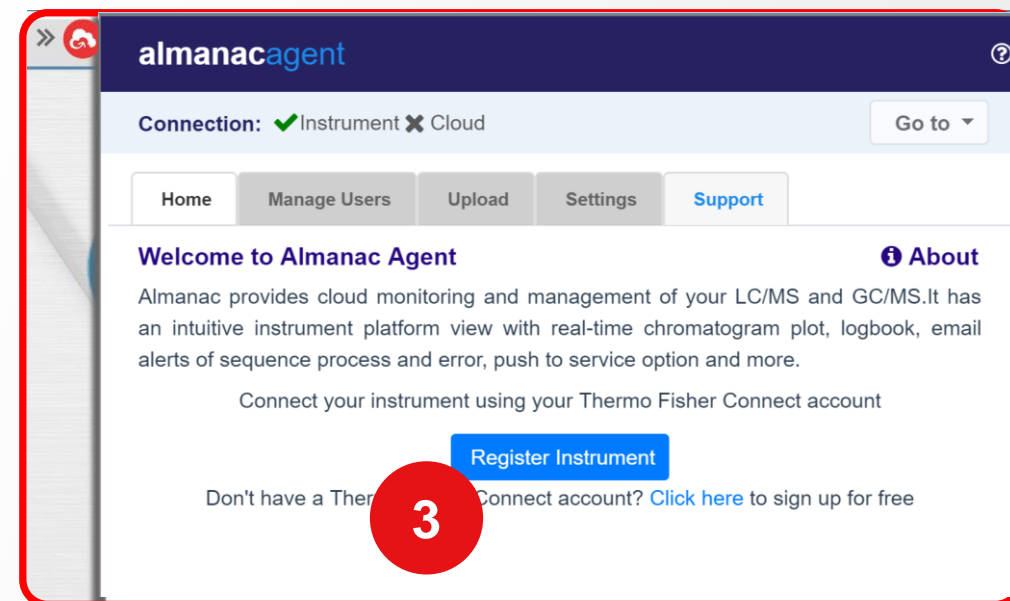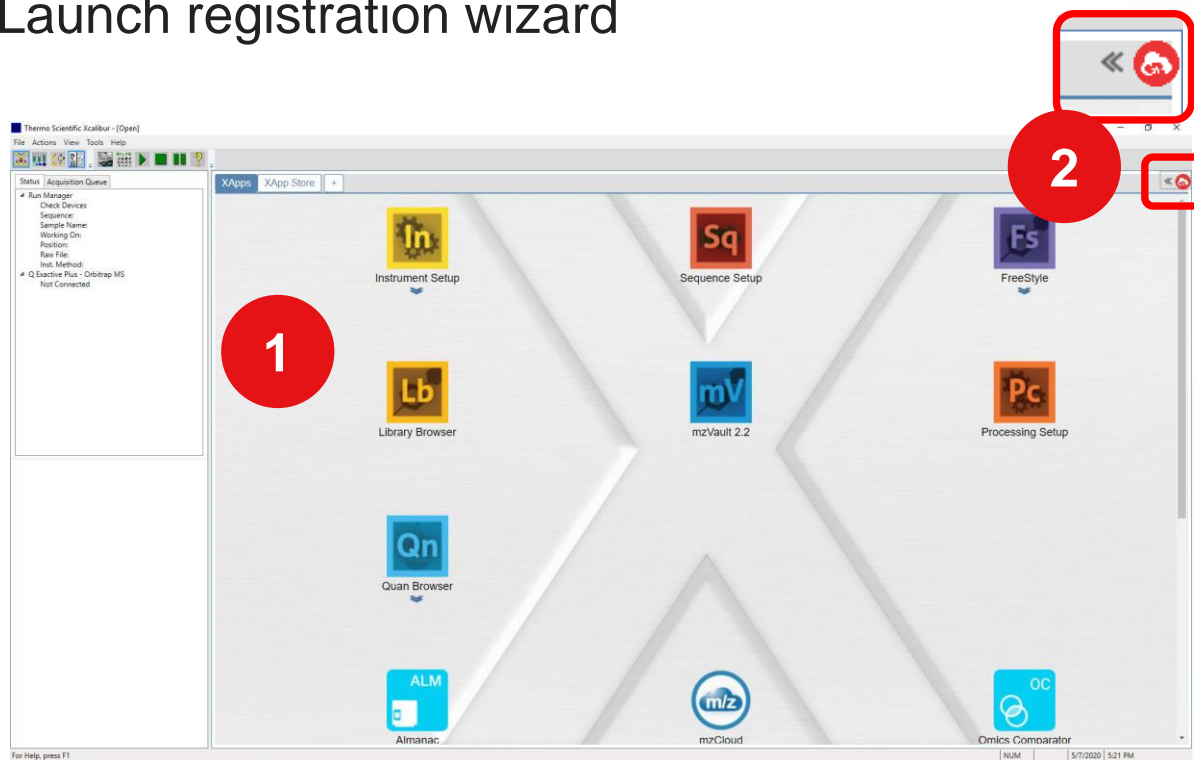

# IAPI Licensing

Via online, click-through agreement

# Step 1: Log into Almanac Agent + Thermo Fisher Cloud Account

## Orbitrap IAPI Online Licensing


- Log into Almanac Agent (see previous section on Almanac for more information)



# Step 2: Open 'About Tune' dialog box in Tune Application

## Orbitrap IAPI Online Licensing

- Settings > About Tune

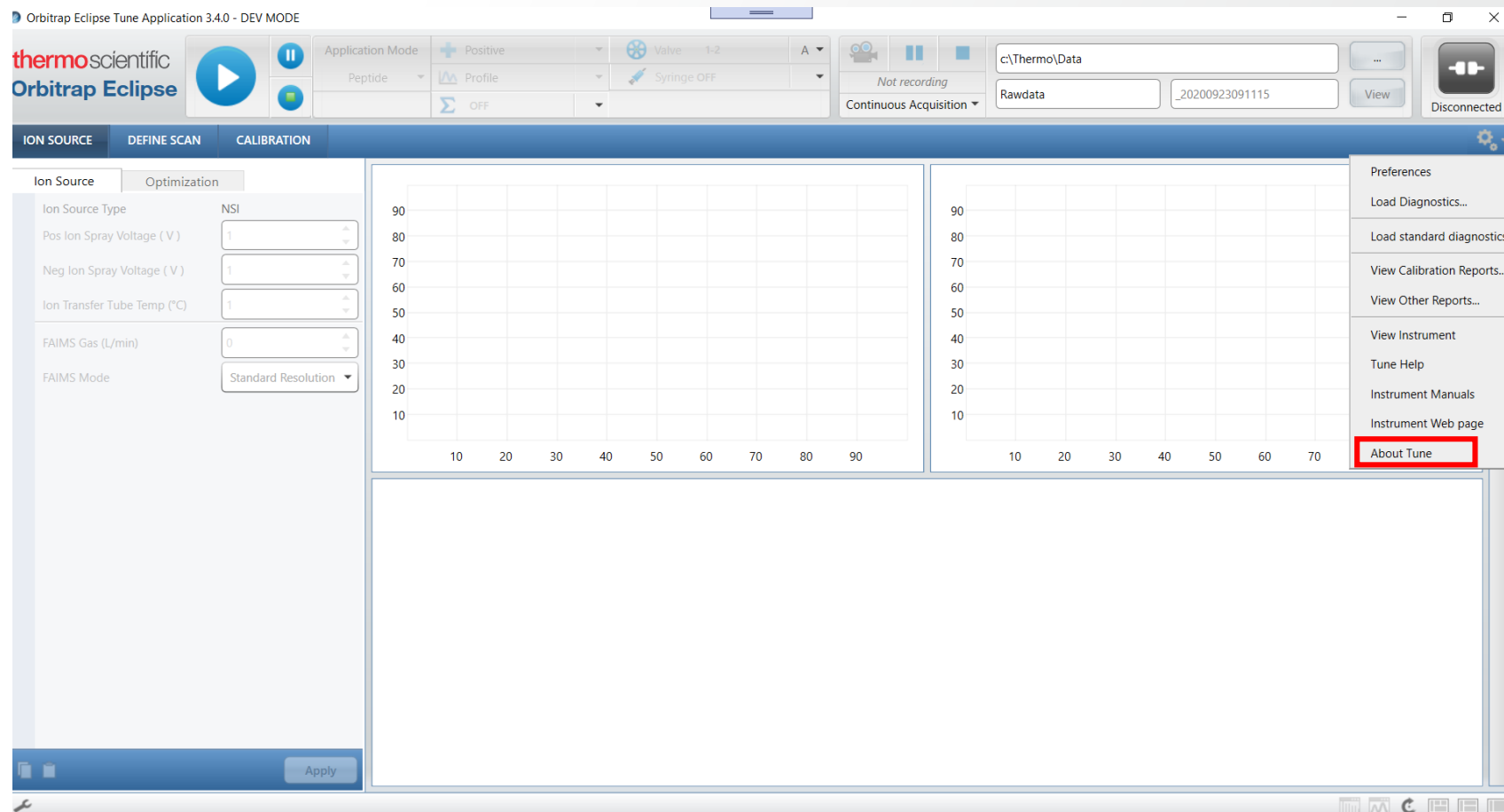




# Step 3: Select Get API License

## Orbitrap IAPI Online Licensing

- About Tune > “License...” dropdown > Get API License

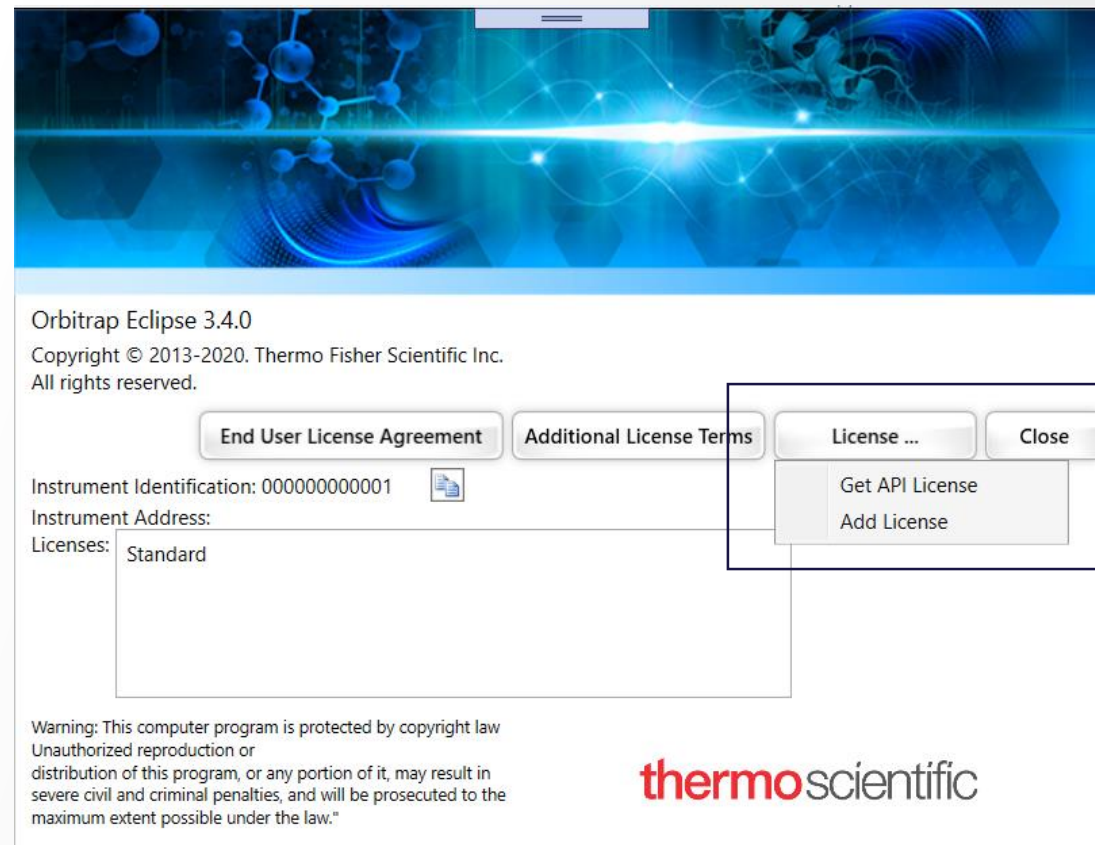

# Step 4: Select IAPI license type

## Orbitrap IAPI Online Licensing

- The dialog shows your instrument information and your account ID
- IAPI license type is the only direct license currently available
- Select “Continue”

The screenshot shows a 'Request License' dialog box with the following content:

Request License

Review the Instrument and User information. Select a License Type to Continue.

Instrument Identification	000000000001
Serial Number	
User Id	yu.zhou2@thermofisher.com
License Type	<div>IAPI</div>

At the bottom of the dialog are three buttons: 'Back', 'Continue', and 'Close'.

# Step 5: Read and Accept License Agreement, Apply Key

## Orbitrap IAPI Online Licensing

- Read through the license agreement and scroll to the bottom to review the instrument details
- Check “I Agree”
- Select “Generate E-Token”
- Check your registered email account for the e-token from [ThermoFisherCloud@thermofisher.com](mailto:ThermoFisherCloud@thermofisher.com)
- Copy and paste the token from the email to the text box area
- Select “Continue”

The screenshot shows a web browser window titled "Request License". The main content area displays the "INSTRUMENT APPLICATION PROGRAMMING INTERFACE (IAPI) LICENSE AGREEMENT". The text of the agreement is visible, including a definition section. A sidebar on the right contains a "Review the License Agreement" section with a checkbox for "I Agree." (callout 2), a scroll bar, and a "Generate E-Token" button (callout 3). Below the sidebar is a text box labeled "Apply E-Token here to continue:" (callout 4). At the bottom of the main content area are "Back", "Continue", and "Close" buttons. Callout 1 points to the "1. DEFINITIONS" section of the agreement.

ThermoFisherCloud@thermofisher.com  
Thermo Fisher Scientific License Registration - Registration E-Token

# Step 6: Activate License & Save Agreement Copy

## Orbitrap IAPI Online Licensing

- Once license is granted, select “Activate License” to apply it
- Select “Save License Agreement” to download the License agreement in PDF format
- Same agreement is also available in your Thermo Fisher account login on <https://apps.thermofisher.com/apps/spa/#/dataconnect>

